


JUGULAR CATHETER SET – FOR EXTENDED USE

WITH PEEL-AWAY® SHEATH INTRODUCTION TECHNIQUE

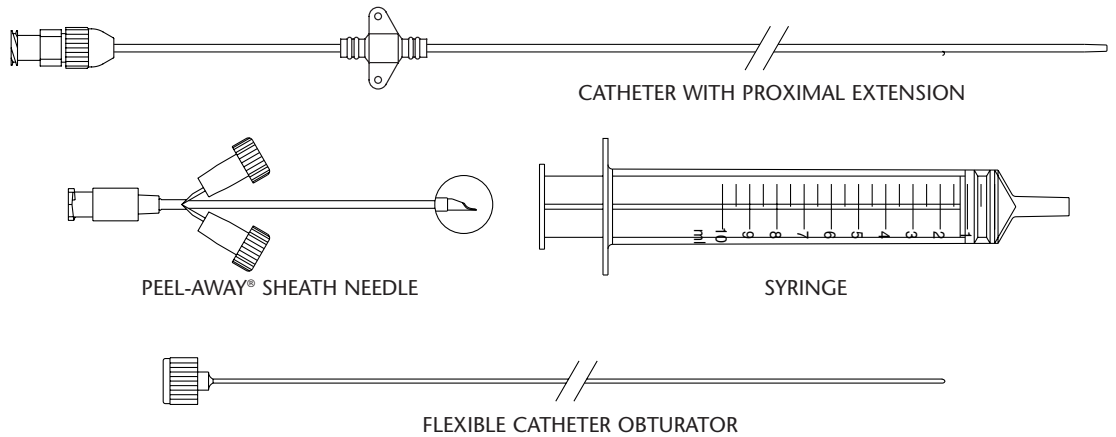

DESCRIPTION

- Central venous catheter introduced using a Peel-Away® sheath needle.
- The Peel-Away® sheath preloaded on the access needle allows percutaneous introduction of the venous catheter, following direct venipuncture and removal of the needle cannula.
- After positioning of the catheter, the sheath is peeled away. Catheter can be trimmed if shorter length is required.
- Proximal extension assists with catheter manipulations after placement.
- Intended for ONE TIME USE.
- Intended for VETERINARY USE ONLY.

ORDER No.

Remarks

V-PPC-5.0U-20-PE5.0

Set consists of a 5 French polyurethane catheter 20 cm (8") long with 5.0 cm proximal catheter extension, a matched flexible obturator, a 15 gauge 7.0 cm long Peel-Away® sheath needle and a 10 ml syringe.

V-PPC-7.0U-25-PE5.0

Set consists of a 7 French polyurethane catheter 25 cm (10") long with 5.0 cm proximal catheter extension, a matched flexible obturator, a 13 gauge 7.0 cm long Peel-Away® sheath needle and a 10 ml syringe.

V-PPC-8.5U-25-PE5.0

Set consists of a 8.5 French polyurethane catheter 25 cm (10") long with 5.0 cm proximal catheter extension, a matched flexible obturator, a 12 gauge 7.0 cm long Peel-Away® sheath needle and a 10 ml syringe.

V-PPC-10.2U-25-PE5.0

Set consists of a 10.2 French polyurethane catheter 25 cm (10") long with 5.0 cm proximal catheter extension, a matched flexible obturator, a 10 gauge 7.0 cm long Peel-Away® sheath needle and a 10 ml syringe.