


OCULAR LAVAGE KIT

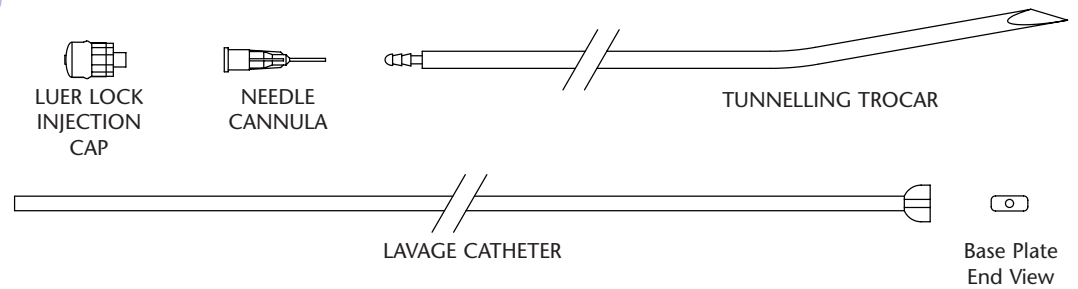

DESCRIPTION

- A silicone endhole catheter used for the controlled irrigation of therapeutic solutions over the eye surface.
- The trocar is used to draw the catheter through the eyelid. The catheter base plate is placed sub-palpebrally and the catheter secured above the exit hole.
- Intended for VETERINARY USE ONLY.

ORDER No.

Remarks

V-POLK5-90

Kit includes 5 French lavage catheter, 90 cm (36") long with 0.76 mm I.D. and 1.7 mm O.D., trocar, blunt tipped needle, needle injection cap and scalpel.

V-POLK8-90

Kit includes 8 French lavage catheter, 90 cm (36") long with 1.5 mm I.D. and 2.5 mm O.D., trocar, blunt tipped needle, needle injection cap and scalpel.

REFERENCES

Veterinary Medicine, May 1992, p460-471
James V Schoster, DVM, Dipl ACVO, University of Wisconsin, Madison, Wisconsin, USA.